

7 COLORS FOR
DIFFERENT TREATMENTS



7 WAVELENGTH. PHOTODYNAMIC THERAPY[®]

Dr. bioled[®]
FACIAL MASK

Distribuido por:




Inlab,S.L. 28430 Alpedrete. Madrid. España. www.inlabmed.es



Dr. bioled[®]
FACIAL MASK

Rechargeable and portable
LED phototherapy device.

Uniform light for facial
cosmetic treatment.

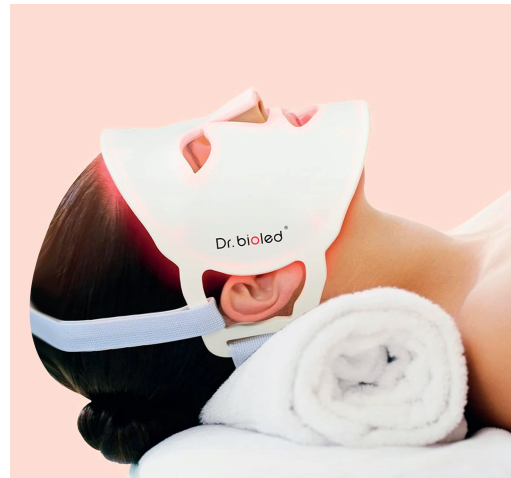

Weight: 400 G.

Dr. bioled®

It is a flexible silicone LED Mask that contains multiple light-emitting diodes. The LED Mask is placed on the face and secured with an adjustable velcro strap.

Dr. bioled® LED mask uses 7 wavelengths of light, which stimulate collagen and elastin producing cells.

EFFECTS

- Treatment that does not produce pain, heat or irritation.
- Improves the texture, appearance and luminosity of the skin.
- Maintenance: 1 weekly session.



Evens skin tone




Helps treat acne



Cellular activator



Smoothing of the skin surface

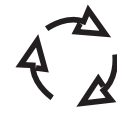

Regulation of oily skin


PRECAUTIONS

- Do not use on pregnant women, children, people with epilepsy, thyroid or dermatological problems, or people allergic to light or medications.
- Securely adjust the strap around your head, before adjusting light intensity, time and color.
- Lying on a bed or sofa with the eyes closed is recommended. Close your eyes and cover them with protective glasses or cotton.
- Avoid eye contact with the light surface or look at it.
- The luminescent film inside the circuit cannot be bent or crushed.
- Do not wet the mask.
- Clean it with a damp cloth. Do not use corrosive chemical solutions
- Avoid sunlight on the mask.
- If you notice any uncomfortable effects, remove the LED mask.



INSTRUCTIONS OF USE

- 1 Place the LED Mask on your face, secure with the Velcro and press power button. Next select the color of the LED light.
- 2 Press the + button to add treatment time. Press - to reduce treatment time. 15-30 minutes treatment time is suggested.
- 3 Press the "sun" icon button for more power or the small sun icon to reduce power.
- 4 After the treatment, apply moisturizing cream and sunscreen.
- 5 Clean the LED Mask with a water-based wipe or damp cloth.

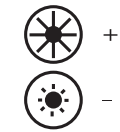

ACTION

7 COLORS FOR DIFFERENT TREATMENTS

- **RED LIGHT: Elasticity**
Stimulates the production of collagen and elastin. Repairs damaged tissue and helps it to heal. Improves the appearance of blemishes and bags under the eyes. Reduces fine expression lines and enlarged pores.
- **BLUE LIGHT: Anti-acne.** Reduces redness. Penetrates deep within the skin. Reduces enlarged pores and helps to prevent the emergence of spots. Generates biostimulation that enhances the natural repair mechanisms, achieving an improvement in skin appearance.
- **GREEN LIGHT: Anti-blemish**
A balance of green light can improve skin blemishes. Reduces the appearance of fine expression lines. Helps to regenerate nutrients in the skin.
- **YELLOW LIGHT: Radiance**
Improves oxygen exchange in the cells, replenishing energy and improving lymphatic drainage. Smooths rough skin and improves wrinkles, making the skin look radiant.
- **PURPLE LIGHT: Healing**
By using red and blue light, two phytotherapeutic effects are combined. Repairs marks and alleviates redness.
- **LIGHT BLUE LIGHT: Unclogging.** Moisturises. It may gradually improve cellular energy.
- **LASER LIGHT (WHITE): Anti-aging**
Penetrates deep within the skin, speeding up skin tissue regeneration. Improves expression lines appearance and helps to firm loose skin.

